


### MORTUARY CHAMBER 6 BODIES - 2 DOORS

Cod. CAM030602

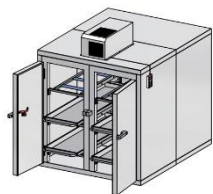

Disclaimer:  
Images are purely  
indicative. The features  
and the accessories of the  
products may differ  
according to model and  
application.

#### Characteristics

<b>Temperature range</b>	+2 / +8 °C <i>available also with -10 -20 °C</i>
<b>Capacity</b>	6 bodies on 3 levels
<b>External Dimensions (WxDxH)</b>	cm 191x226x250 <i>Dimensions without monoblock: cm. 191x226x211; weight 850 kg</i>
<b>Structure and Insulation</b>	White coated steel inside and outside; insulation polyurethane thickness 80 mm
<b>Door</b>	n.2 doors
<b>Key Lock</b>	
<b>Light</b>	
<b>Internal equipment</b>	n.6 steel body trays
<b>Controller</b>	Digital electronic thermostat
<b>Alarms</b>	Acoustic and Visual for temperature deviations and failures
<b>Refrigeration Type</b>	Ventilated
<b>Defrost</b>	Automatic with automatic evaporation of the condensate water

#### Technical data

<b>Details</b>	<ul style="list-style-type: none"> <li>• Refrigerant type: R452A, CFC free</li> <li>• Mains Voltage 220V/1/50Hz (-20°C model: 380V/3/50Hz)</li> <li>• HP 1,2 E ;kW 0,7 (-20°C model: HP 2 E; kW 1,5)</li> <li>• Refrigerating capacity: 2245 W (-20°C model: 1686 W)</li> <li>• Ambient Temperature: +10/+32 °C; Humidity 60%</li> <li>• Connector: Schuko</li> </ul>
<b>In accordance with</b>	2006/42/CE, 2014/30/CE, 2014/35/CE and following modifications
<b>Package</b>	Chamber to be installed on site, packed in cellophane on pallet, monoblock in strong cardboard on wooden pallet
<b>Note</b>	Available also with a negative temperature range (-10° and -20°C)

Internal covering in stainless steel

External covering in stainless steel

Small hatches, instead of the big door

### DMLP UNIT

External module with DMLP DIGITAL MONITOR:  
Monitor/alarm/recording system with rechargeable battery  
- Temperatures and alarms recorder on removable Secure Digital (SD Back-Up memory): it can store 5 years data, recording temperature and alarms every minute.  
The SD card can be connected to any PC: a software (Capture) is included, so that you can download, save and print all data from your PC  
- Battery backup: ensures alarms and recording of temperatures even in case of power failure. The battery system lasts about 12 hours.  
- Dry contact for remote connection + RS485 port

### DMLP PRINTER

Additional body tray

Hydraulic trolley cm 30-165

Hydraulic trolley with 12V electrical adjustment cm 30-165

Cover for hidraulic trolley

60 Hz Frequency

### FIOCCHETTI SCIENTIFIC S.R.L.

Headquarter: Via Panagulis, 48 - 42045 Luzzara (RE) Italy  
VAT No. IT00934960352

Phone +39 0522 976232 - Fax +39 0522 976028

[www.fiocchetti.com](http://www.fiocchetti.com) - Email: [info@fiocchetti.it](mailto:info@fiocchetti.it) - Pec: [fiocchetti@pec.it](mailto:fiocchetti@pec.it)